

#### KEY NOTES

1. MOUNT A5145, HOOK, AT 1372 [4'-6"] AFF.
2. INCREASE THIS DIMENSION AS NEEDED TO COMPLY WITH ADA IF THE DOOR HAS A CLOSER OR IF 610 [2'-0"] CLEARANCE IS NOT AVAILABLE AT THE STRIKE SIDE OF THE DOOR.
3. COORDINATE LOCATION OF M1620, PATIENT CHART HOLDER, AND ROOM SIGNAGE.
4. LENGTH OF E1500, WALL MOUNTED RAIL, TO BE DETERMINED BY MANUFACTURER/INSTALLER.
5. MOUNT A5145, HOOK, AT 1219 [4'-0"] AFF. PROVIDE 2-A5145'S FOR NAVY PROJECTS ONLY. MOUNT 1-A5145 AT 1219 [4'-0"] AFF AND 1-A5145 AT 1524 [5'-0"] AFF.
6. P3200, ELECTRONIC LAVATORY, MAY BE USED INSTEAD OF P3100, LAVATORY. FOR THE P3200 PROVIDE EITHER BATTERY POWER OR HARD WIRED CONNECTION.

## OUTPATIENT DERMATOLOGY, CRYOTHERAPY (OPDC1)

**3.15.01**

SPECIALTY MEDICAL CLINICS


FLOOR/EQUIPMENT PLAN 11.15 NSM (120 NSF)

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

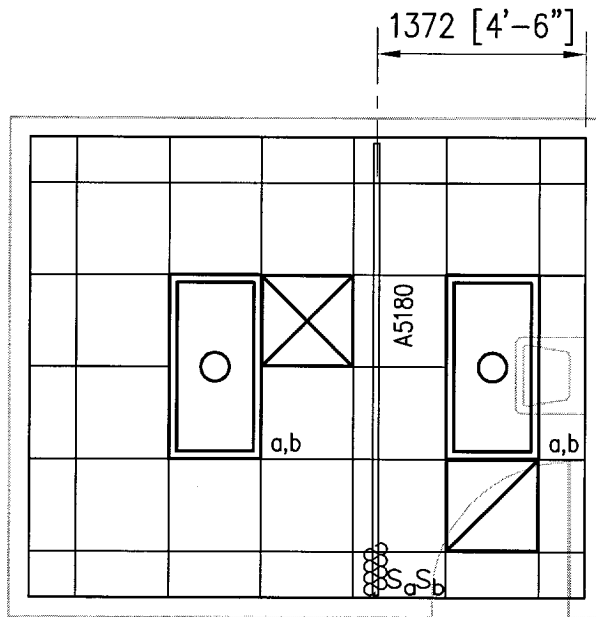

# **OUTPATIENT DERMATOLOGY, CRYOTHERAPY (OPDC1)**

**3.15.01**

SPECIALTY MEDICAL CLINICS


CEILING PLAN

DATE: 01/28/2002


0 1000 2000




1:50 SCALE (MILLIMETERS)

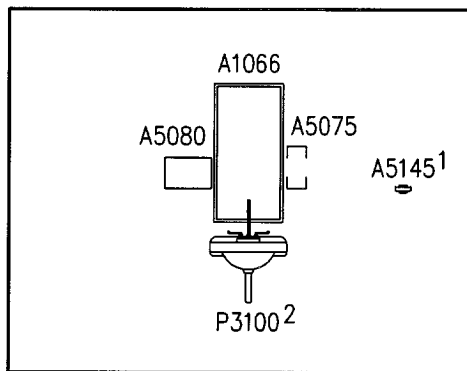

**EAST**

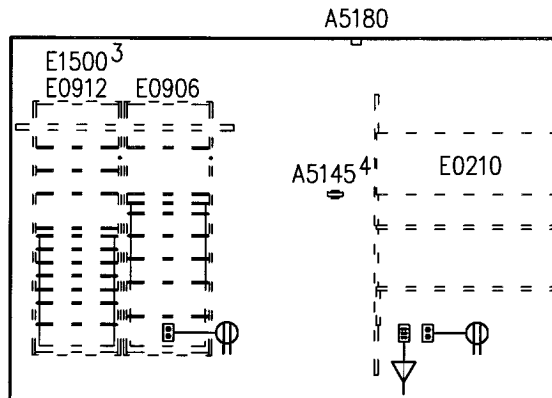

**NORTH**

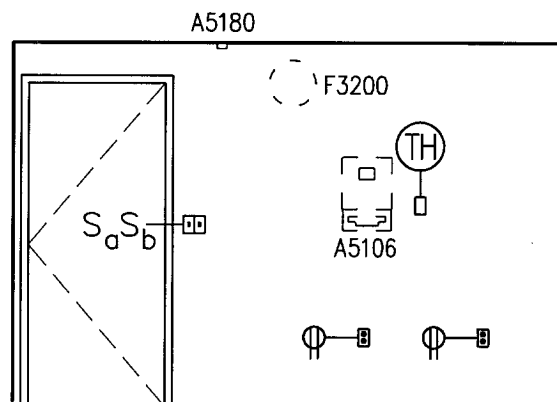

**SOUTH**

**NOTES**

1. SEE TYPICAL EQUIPMENT MOUNTING HEIGHTS BEGINNING ON SHEET 7.01.01.

**KEY NOTES**

1. MOUNT A5145, HOOK, AT 1219 [4'-0"] AFF. PROVIDE 2-A5145'S FOR NAVY PROJECTS ONLY. MOUNT 1-A5145 AT 1219 [4'-0"] AFF AND 1-A5145 AT 1524 [5'-0"] AFF.
2. P3200, ELECTRONIC LAVATORY, MAY BE USED INSTEAD OF P3100, LAVATORY. FOR THE P3200 PROVIDE EITHER BATTERY POWER OR HARD WIRED CONNECTION.
3. LENGTH OF E1500, WALL MOUNTED RAIL, TO BE DETERMINED BY MANUFACTURER/INSTALLER.
4. MOUNT A5145, HOOK, AT 1372 [4'-6"] AFF.

# **OUTPATIENT DERMATOLOGY, CRYOTHERAPY (OPDC1)**

**3.15.01**

SPECIALTY MEDICAL CLINICS

ELEVATIONS

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

## Medical Equipment Room Guide Plates


### 3.15 Specialty Medical Clinics

**Plate Number** 3.15.01

**Room Name** OUTPATIENT DERMATOLOGY, CRYOTHERAPY

**SEPS #** OPDC1

JSN	Nomenclature	Qty	Army	Navy	AF	Utilities					
			Cat	Cat	Cat	1	2	3	4	5	6
A1066	Mirror, Float Glass, With SS Frame, 36X18	1	A	A	A	.	.	.	.	.	J
A5075	Dispenser, Soap, Disposable	1	C	C	C	.	.	.	.	.	.
A5080	Dispenser, Paper Towel, SS, Surface Mounted	1	A	A	A	.	.	.	.	.	J
A5106	Waste Disposal Unit, Sharps w/Glove Dispenser	1	C	A	C	.	.	.	.	.	.
A5145	Hook, Garment, Double, SS, Surface Mounted	2	A	A	A	.	.	.	.	.	J
A5180	Track, Cubicle, Surface Mounted, With Curtain	1	A	A	A	.	.	.	.	.	J
E0210	Worksurface, w/Overhead Cab, Wall Mtd, 48" W	1	C	C	C	.	A	.	.	.	J
E0906	Locker, Supply, General, Wall Mtd, 23"W x 20"D	1	C	C	C	.	.	.	.	.	.
E0912	Locker, Supply, Med Surg, Wall Mtd, 23"W x 20"D	1	C	C	C	.	.	.	.	.	.
E0948	Cart, General Storage, Mobile, 42"H x 32"W x 22"D	1	C	C	C	.	.	.	.	.	.
E1500	Rail, MOD, W/MNTD, HX144XD	1	C	C	C	.	.	.	.	.	J
F0205	Chair, Side With Arms	1	C	C	C	.	.	.	.	.	.
F0280	Chair, Swivel, Low Back	1	C	C	C	.	.	.	.	.	.
F0340	Stool, Self Adjusting	1	C	C	C	.	.	.	.	.	.
F2000	Basket, Wastepaper, Round, Metal, 18 H x 16 Dia.	1	C	C	C	.	.	.	.	.	.
F3200	Clock, Battery, 12" Diameter	1	C	C	C	.	.	.	.	.	.
M1620	Holder, Chart, Patient, Wall or Door Mounted	1	A	A	A	.	.	.	.	.	.
M1800	Computer, Microprocessing	1	C	C	C	.	A	.	.	.	.
M3072	Frame, Infectious Waste Bag w/Lid	1	C	C	C	.	.	.	.	.	.
M4925	Chair, Exam/Treatment, With Motor	1	C	C	C	.	A	.	.	.	.
M7410	Exam Light, Ultraviolet (Woods)	1	C	C	C	.	A	.	.	.	.
M8810	Stand, Mayo	1	C	C	C	.	.	.	.	.	.
P3100	Lavatory, Vitreous China, Slab Type	1	A	A	A	D	.	.	.	.	.

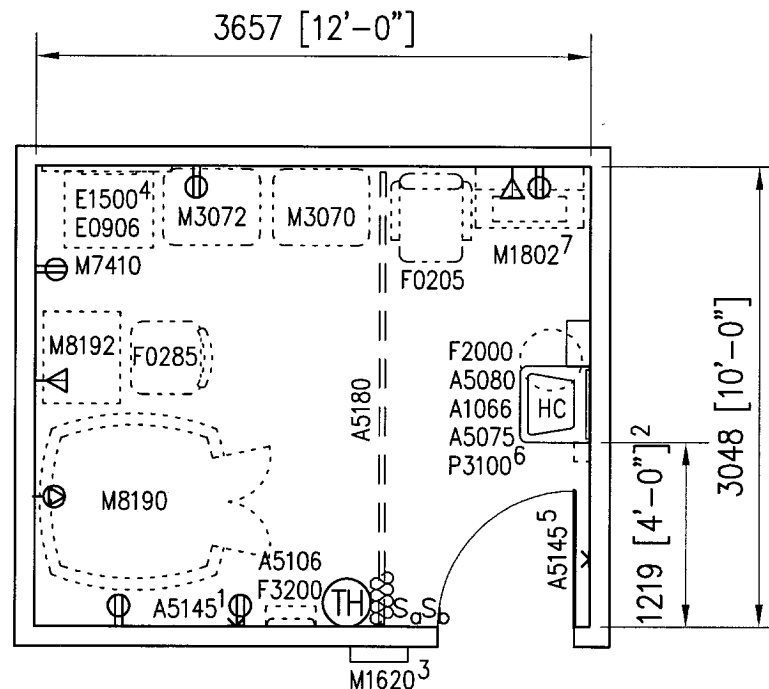

#### KEY NOTES

1. MOUNT A5145, HOOK, AT 1372 [4'-6"] AFF.
2. INCREASE THIS DIMENSION AS NEEDED TO COMPLY WITH ADA IF THE DOOR HAS A CLOSER OR IF 610 [2'-0"] CLEARANCE IS NOT AVAILABLE AT THE STRIKE SIDE OF THE DOOR.
3. COORDINATE LOCATION OF M1620, PATIENT CHART HOLDER, AND ROOM SIGNAGE.
4. LENGTH OF E1500, WALL MOUNTED RAIL, TO BE DETERMINED BY MANUFACTURER/INSTALLER.
5. MOUNT A5145, HOOK, AT 1219 [4'-0"] AFF. PROVIDE 2-A5145'S FOR NAVY PROJECTS ONLY. MOUNT 1-A5145 AT 1219 [4'-0"] AFF AND 1-A5145 AT 1524 [5'-0"] AFF.
6. P3200, ELECTRONIC LAVATORY, MAY BE USED INSTEAD OF P3100, LAVATORY. FOR THE P3200 PROVIDE EITHER BATTERY POWER OR HARD WIRED CONNECTION.
7. COORDINATE LOCATIONS OF ELECTRICAL AND DATA OUTLETS WITH M1802, WALL MOUNTED, FOLD DOWN COMPUTER.

## OUTPATIENT DERMATOLOGY, UV BOOTH (OPDU1) 3.15.02

SPECIALTY MEDICAL CLINICS


FLOOR/EQUIPMENT PLAN 11.15 NSM (120 NSF)

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

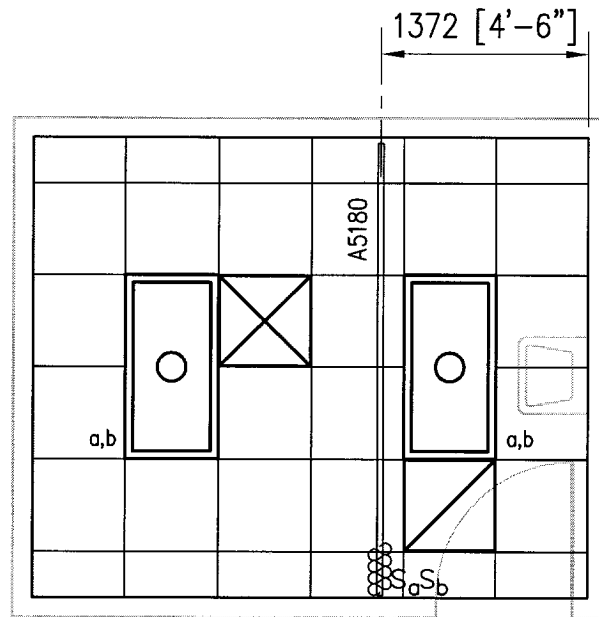

# **OUTPATIENT DERMATOLOGY, UV BOOTH (OPDU1) 3.15.02**

SPECIALTY MEDICAL CLINICS


CEILING PLAN

DATE: 01/28/2002

0 1000 2000




1:50 SCALE (MILLIMETERS)

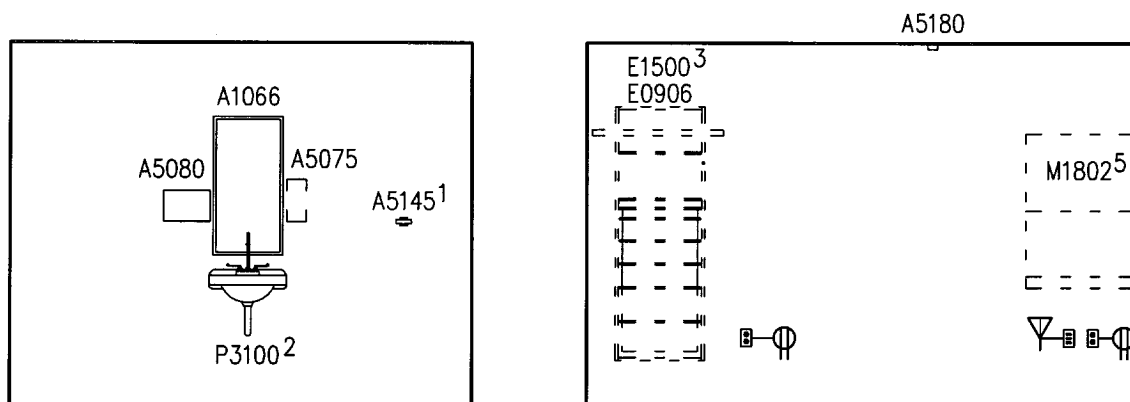

**EAST**

**NORTH**

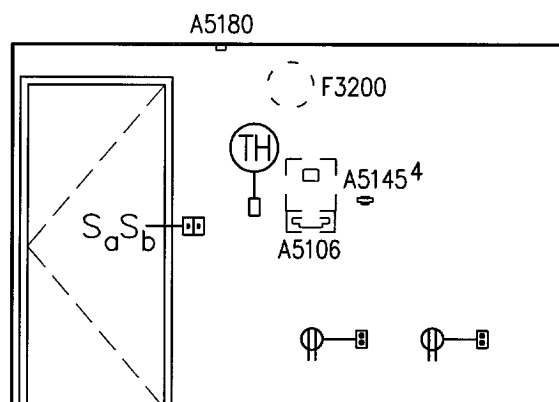

**SOUTH**

**NOTES**

1. SEE TYPICAL EQUIPMENT MOUNTING HEIGHTS BEGINNING ON SHEET 7.01.01.

**KEY NOTES**

1. MOUNT A5145, HOOK, AT 1219 [4'-0"] AFF. PROVIDE 2-A5145'S FOR NAVY PROJECTS ONLY. MOUNT 1-A5145 AT 1219 [4'-0"] AFF AND 1-A5145 AT 1524 [5'-0"] AFF.
2. P3200, ELECTRONIC LAVATORY, MAY BE USED INSTEAD OF P3100, LAVATORY. FOR THE P3200 PROVIDE EITHER BATTERY POWER OR HARD WIRED CONNECTION.
3. LENGTH OF E1500, WALL MOUNTED RAIL, TO BE DETERMINED BY MANUFACTURER/INSTALLER.
4. MOUNT A5145, HOOK, AT 1372 [4'-6"] AFF.
5. COORDINATE LOCATIONS OF ELECTRICAL AND DATA OUTLETS WITH M1802, WALL MOUNTED, FOLD DOWN COMPUTER.

# **OUTPATIENT DERMATOLOGY, UV BOOTH (OPDU1) 3.15.02**

SPECIALTY MEDICAL CLINICS

ELEVATIONS

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

## Medical Equipment Room Guide Plates

### 3.15 Specialty Medical Clinics


**Plate Number** 3.15.02

**Room Name** OUTPATIENT DERMATOLOGY, UV BOOTH

**SEPS #** OPDU1

JSN	Nomenclature	Qty	Army	Navy	AF	Utilities					
			Cat	Cat	Cat	1	2	3	4	5	6
A1066	Mirror, Float Glass, With SS Frame, 36X18	1	A	A	A	.	.	.	.	.	J
A5075	Dispenser, Soap, Disposable	1	C	C	C	.	.	.	.	.	.
A5080	Dispenser, Paper Towel, SS, Surface Mounted	1	A	A	A	.	.	.	.	.	J
A5106	Waste Disposal Unit, Sharps w/Glove Dispenser	1	C	A	C	.	.	.	.	.	.
A5145	Hook, Garment, Double, SS, Surface Mounted	2	A	A	A	.	.	.	.	.	J
A5180	Track, Cubicle, Surface Mounted, With Curtain	1	A	A	A	.	.	.	.	.	J
E0906	Locker, Supply, General, Wall Mtd, 23"W x 20"D	1	C	C	C	.	.	.	.	.	.
E1500	Rail, MOD, W/MNTD, HX144XD	1	C	C	C	.	.	.	.	.	J
F0205	Chair, Side With Arms	1	C	C	C	.	.	.	.	.	.
F0285	Chair, Secretarial, Tilt Back, Adjustable Height	1	C	C	C	.	.	.	.	.	.
F2000	Basket, Wastepaper, Round, Metal, 18 H x 16 Dia.	1	C	C	C	.	.	.	.	.	.
F3200	Clock, Battery, 12" Diameter	1	C	C	C	.	.	.	.	.	.
M1620	Holder, Chart, Patient, Wall or Door Mounted	1	A	A	A	.	.	.	.	.	.
M1802	Wall Mounted, Fold Down Computer	1	C	C	C	.	A	.	.	.	J
M3070	Hamper, Linen, Mobile, w/Lid	1	C	C	C	.	.	.	.	.	.
M3072	Frame, Infectious Waste Bag w/Lid	1	C	C	C	.	.	.	.	.	.
M7410	Exam Light, Ultraviolet (Woods)	1	C	C	C	.	A	.	.	.	.
M8190	Booth, Treatment, Ultraviolet	1	C	C	C	.	C	.	.	.	A
M8192	Booth, Treatment, Ultraviolet, Hand-Foot	1	C	C	C	.	A	.	.	.	.
P3100	Lavatory, Vitreous China, Slab Type	1	A	A	A	D	.	.	.	.	.



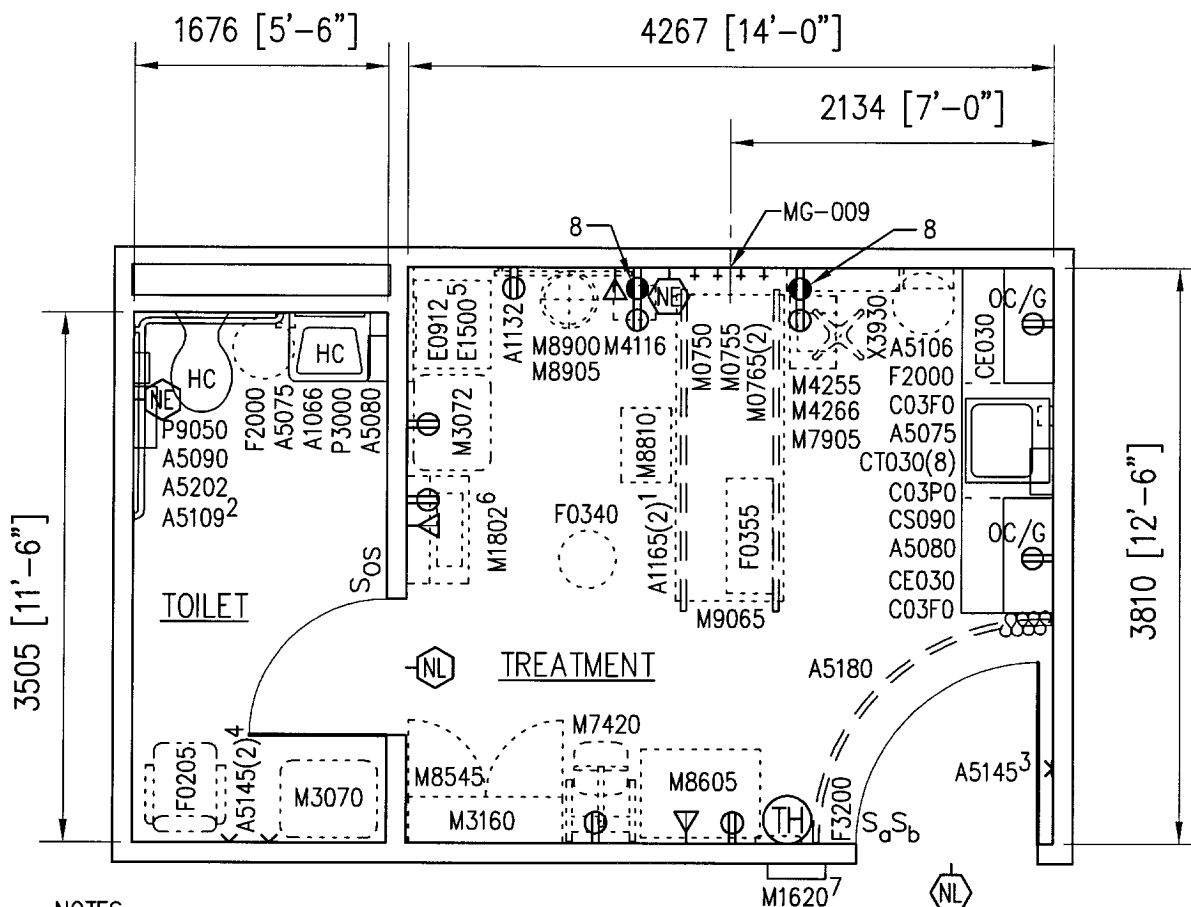

#### NOTES

1. PROVIDE ONE EACH E0954, EMERGENCY CART, AND M7660, DEFIBRILLATOR, PER CLINIC.

#### KEY NOTES

1. PROVIDE A1165, IV TRACKS, FOR NAVY PROJECTS ONLY.
2. PROVIDE A5109, GRAB BAR, THAT COMPLIES WITH ADA.
3. MOUNT A5145, HOOK, AT 1372 [4'-6"] AFF. PROVIDE 2-A5145'S FOR NAVY PROJECTS ONLY. MOUNT 1-A5145 AT 1372 [4'-6"] AFF AND 1-A5145 AT 1524 [5'-0"] AFF.
4. MOUNT 1-A5145 AT 1372 [4'-6"] AFF AND 1-A5145 AT 1524 [5'-0"] AFF.
5. LENGTH OF E1500, WALL MOUNTED RAIL, TO BE DETERMINED BY MANUFACTURER/INSTALLER.
6. COORDINATE LOCATIONS OF ELECTRICAL AND DATA OUTLETS WITH M1802, WALL MOUNTED, FOLD DOWN COMPUTER.
7. COORDINATE LOCATION OF M1620, PATIENT CHART HOLDER, AND ROOM SIGNAGE.
8. UTILITIES AND DEVICES INDICATED ARE FOR A CLINIC OR HOSPITAL WITH AN ALTERNATE POWER SOURCE.

## TREATMENT ROOM, PROCTOSCOPIC EXAM (TRPE1)

**3.15.03**

SPECIALTY MEDICAL CLINICS

FLOOR/EQUIPMENT PLAN


22.13 NSM (238 NSF)

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

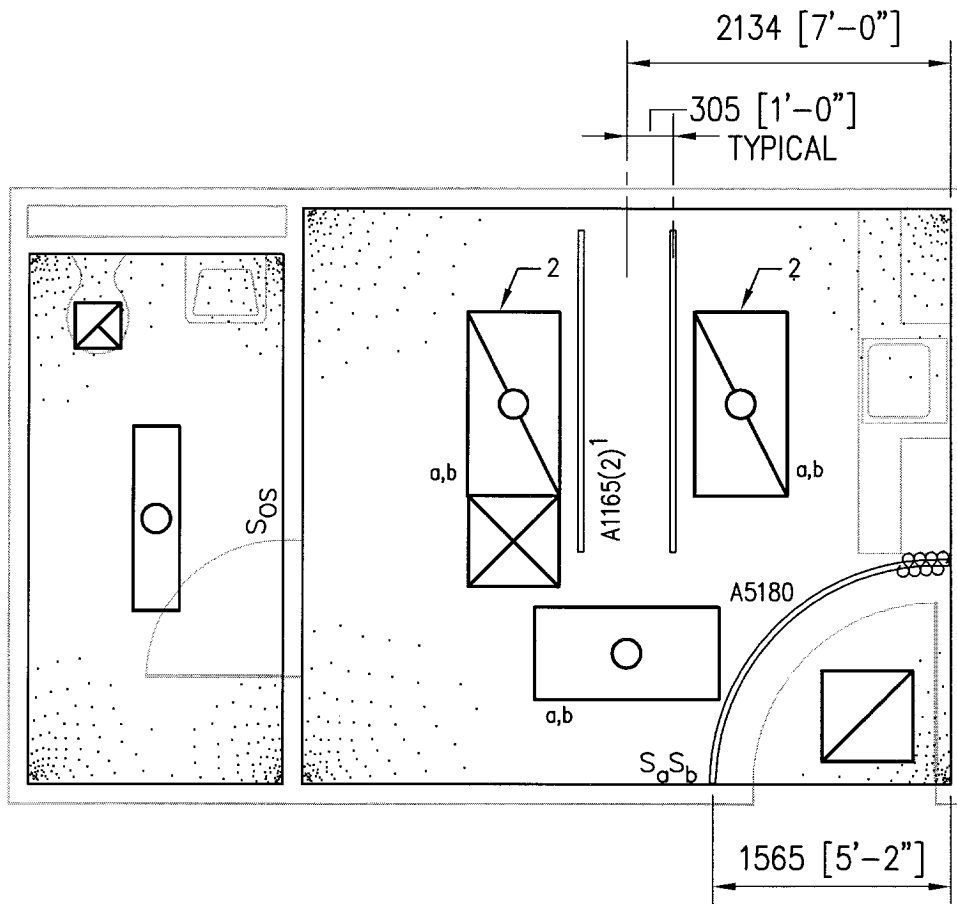

#### KEY NOTES

1. PROVIDE A1165, IV TRACKS, FOR NAVY PROJECTS ONLY.
2. UTILITIES AND DEVICES INDICATED ARE FOR A CLINIC OR HOSPITAL WITH AN ALTERNATE POWER SOURCE. WHERE NO ALTERNATE POWER SOURCE IS AVAILABLE PROVIDE BATTERY BACKUP BALLAST.

## TREATMENT ROOM, PROCTOSCOPIC EXAM (TRPE1)

**3.15.03**

SPECIALTY MEDICAL CLINICS


CEILING PLAN

DATE: 01/28/2002


0 1000 2000




1:50 SCALE (MILLIMETERS)

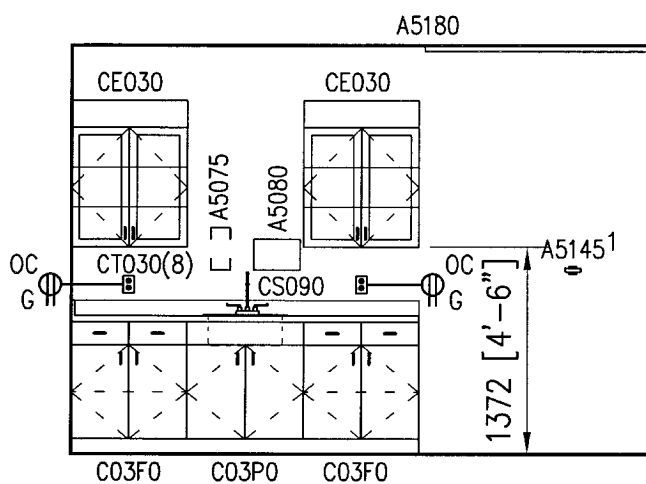

**EAST**

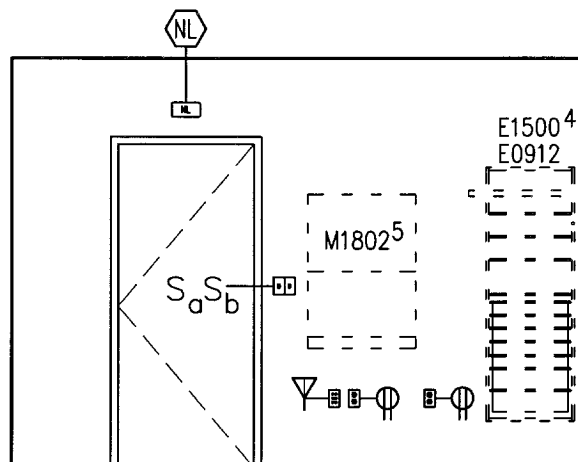

**WEST**

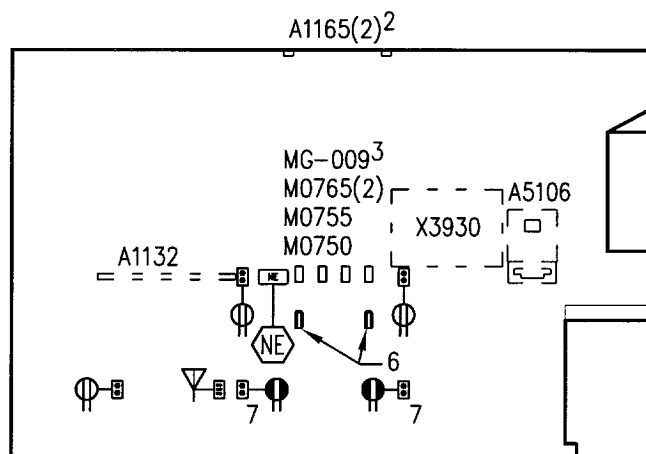

**NORTH**

**NOTES**

1. SEE TYPICAL EQUIPMENT MOUNTING HEIGHTS BEGINNING ON SHEET 7.01.01.

**KEY NOTES**

1. MOUNT A5145, HOOK, AT 1372 [4'-6"] AFF. PROVIDE 2-A5145'S FOR NAVY PROJECTS ONLY. MOUNT 1-A5145 AT 1372 [4'-6"] AFF AND 1-A5145 AT 1524 [5'-0"] AFF.
2. PROVIDE A1165, IV TRACKS, FOR NAVY PROJECTS ONLY.
3. PROVIDE BOTTLED GAS IF BUILDING PIPED MCA AND OX ARE NOT AVAILABLE.
4. LENGTH OF E1500, WALL MOUNTED RAIL, TO BE DETERMINED BY MANUFACTURER/INSTALLER.
5. COORDINATE LOCATIONS OF ELECTRICAL AND DATA OUTLETS WITH M1802, WALL MOUNTED, FOLD DOWN COMPUTER.
6. VACUUM SLIDE.
7. UTILITIES AND DEVICES INDICATED ARE FOR A CLINIC OR HOSPITAL WITH AN ALTERNATE POWER SOURCE.

**TREATMENT ROOM,  
PROCTOSCOPIC EXAM (TRPE1)**

**3.15.03**

SPECIALTY MEDICAL CLINICS


ELEVATIONS

DATE: 01/28/2002


0 1000 2000



1:50 SCALE (MILLIMETERS)

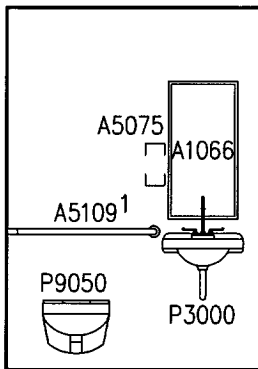

**NORTH** (TOILET)

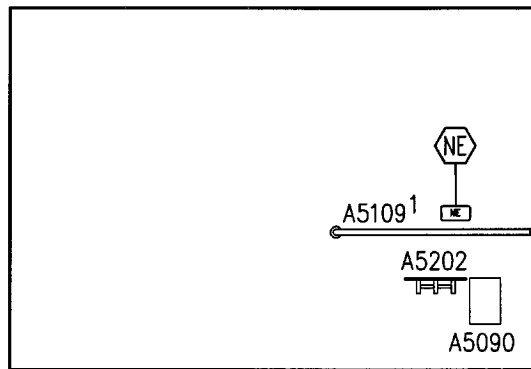

**WEST** (TOILET)

**NOTES**

1. SEE TYPICAL EQUIPMENT MOUNTING HEIGHTS BEGINNING ON SHEET 7.01.01.

**KEY NOTES**

1. PROVIDE A5109, GRAB BAR, THAT COMPLIES WITH ADA.

**TREATMENT ROOM,  
PROCTOSCOPIC EXAM (TRPE1)**

**3.15.03**

SPECIALTY MEDICAL CLINICS

ELEVATIONS

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

# Medical Equipment Room Guide Plates

## 3.15 Specialty Medical Clinics

**Plate Number** 3.15.03

**Room Name** TREATMENT ROOM, PROCTOSCOPIC EXAM (TREATMENT)

**SEPS #** TRPE1

JSN	Nomenclature	Qty	Army	Navy	AF	Utilities					
			Cat	Cat	Cat	1	2	3	4	5	6
A1132	Rail, Accessory Mounting, Length As Required	1	C	C	C	.	.	.	.	.	.
A1165	Track, IV, Ceiling Mounted, 7 Foot	2	A	A	A	.	.	.	.	.	J
A5075	Dispenser, Soap, Disposable	1	C	C	C	.	.	.	.	.	.
A5080	Dispenser, Paper Towel, SS, Surface Mounted	1	A	A	A	.	.	.	.	.	J
A5106	Waste Disposal Unit, Sharps w/Glove Dispenser	1	C	A	C	.	.	.	.	.	.
A5145	Hook, Garment, Double, SS, Surface Mounted	1	A	A	A	.	.	.	.	.	J
A5180	Track, Cubicle, Surface Mounted, With Curtain	1	A	A	A	.	.	.	.	.	J
C03F0	Cabinet, U/C/B, 1 Shelf, 2 Half DR, 2 DO, 36x30x22	2	A	A	A	.	.	.	.	.	.
C03P0	Cabinet, Sink, U/C/B, 2 Door, 36x30x22	1	A	A	A	.	.	.	.	.	.
CE030	Cabinet, W/H, 2 SH, 2 GDO, Sloping Top, 38x30x13	2	A	A	A	.	.	.	.	.	J
CS090	Sink, SS, Single Compartment, 7.5x19x16 ID	1	A	A	A	D	.	.	.	.	.
CT030	Countertop, High Pressure Laminate	8	A	A	A	.	.	.	.	.	.
E0912	Locker, Supply, Med Surg, Wall Mtd, 23"W x 20"D	1	C	C	C	.	.	.	.	.	.
E1500	Rail, MOD, W/MNTD, HX144XD	1	C	C	C	.	.	.	.	.	J
F0205	Chair, Side With Arms	1	C	C	C	.	.	.	.	.	.
F0340	Stool, Self Adjusting	1	C	C	C	.	.	.	.	.	.
F0355	Footstool, Straight	1	C	C	C	.	.	.	.	.	.
F2000	Basket, Wastepaper, Round, Metal, 18 H x 16 Dia.	1	C	C	C	.	.	.	.	.	.
F3200	Clock, Battery, 12" Diameter	1	C	C	C	.	.	.	.	.	.
M0750	Flowmeter, Air, Connect w/50 PSI Supply	1	C	C	C	.	.	K	.	.	.
M0755	Flowmeter, Oxygen, Low Flow	1	C	C	C	.	.	A	.	.	.
M0765	Regulator, Vacuum	2	C	C	C	.	.	B	.	.	.
M1620	Holder, Chart, Patient, Wall or Door Mounted	1	A	A	A	.	.	.	.	.	.
M1802	Wall Mounted, Fold Down Computer	1	C	C	C	.	A	.	.	.	J
M3072	Frame, Infectious Waste Bag w/Lid	1	C	C	C	.	.	.	.	.	.
M3160	Cabinet, Storage, Hanging Scopes	1	C	C	C	.	.	.	.	.	.
M4116	Monitor, Vital Signs	1	C	C	C	.	A	.	.	.	.
M4255	Stand, IV, Adjustable	1	C	C	C	.	.	.	.	.	.
M4266	Pump, Volumetric, Infusion, Multiple Line	1	C	C	C	.	A	.	.	.	.
M7420	Light, Exam, Mobile	1	C	C	C	.	A	.	.	.	.
M7905	Oximeter, Pulse	1	C	C	C	.	A	.	.	.	.
M8545	Sigmoidoscope, Fiberoptic, Therapeutic	1	C	C	C	.	.	.	.	.	.
M8605	Endoscopy Cart, Fiberoptic, Basic	1	C	C	C	.	A	.	.	.	.
M8810	Stand, Mayo	1	C	C	C	.	.	.	.	.	.
M8900	Carriage, Pail, CRS, Without Pail	1	C	C	C	.	.	.	.	.	.
M8905	Pail, Utility, CRS, Without Carriage	1	C	C	C	.	.	.	.	.	.
M9065	Table, Examination, Proctology, Electro/Hydraulic	1	C	C	C	.	A	.	.	.	.
X3930	Illuminator, Film, Double, Wall Mounted, 20x29x6	1	C	C	C	.	A	.	.	.	J

## Medical Equipment Room Guide Plates


### 3.15 Specialty Medical Clinics

**Plate Number** 3.15.03A

**Room Name** TREATMENT ROOM, PROCTOSCOPIC EXAM (TOILET)

**SEPS #** TRPE1

JSN	Nomenclature	Qty	Army	Navy	AF	Utilities					
			Cat	Cat	Cat	1	2	3	4	5	6
A1066	Mirror, Float Glass, With SS Frame, 36X18	1	A	A	A	.	.	.	.	.	J
A5075	Dispenser, Soap, Disposable	1	C	C	C	.	.	.	.	.	.
A5080	Dispenser, Paper Towel, SS, Surface Mounted	1	A	A	A	.	.	.	.	.	J
A5090	Disposal, Sanitary Napkin, SS, Surface Mounted	1	A	A	A	.	.	.	.	.	J
A5109	Grab Bar, 1-1/4" Dia., SS, 2 Wall, W/C Accessible	1	A	A	A	.	.	.	.	.	J
A5145	Hook, Garment, Double, SS, Surface Mounted	2	A	A	A	.	.	.	.	.	J
A5202	Dispenser, Toilet Paper w/Utility Shelf, SS, 2-Roll	1	A	A	A	.	.	.	.	.	J
F2000	Basket, Wastepaper, Round, Metal, 18 H x 16 Dia.	1	C	C	C	.	.	.	.	.	.
M3070	Hamper, Linen, Mobile, w/Lid	1	C	C	C	.	.	.	.	.	.
P3000	Lavatory, Vitreous China, Straight Back	1	A	A	A	D	.	.	.	.	.
P9050	Toilet, Wall Hung, Siphon Jet	1	A	A	A	B	.	.	.	.	.

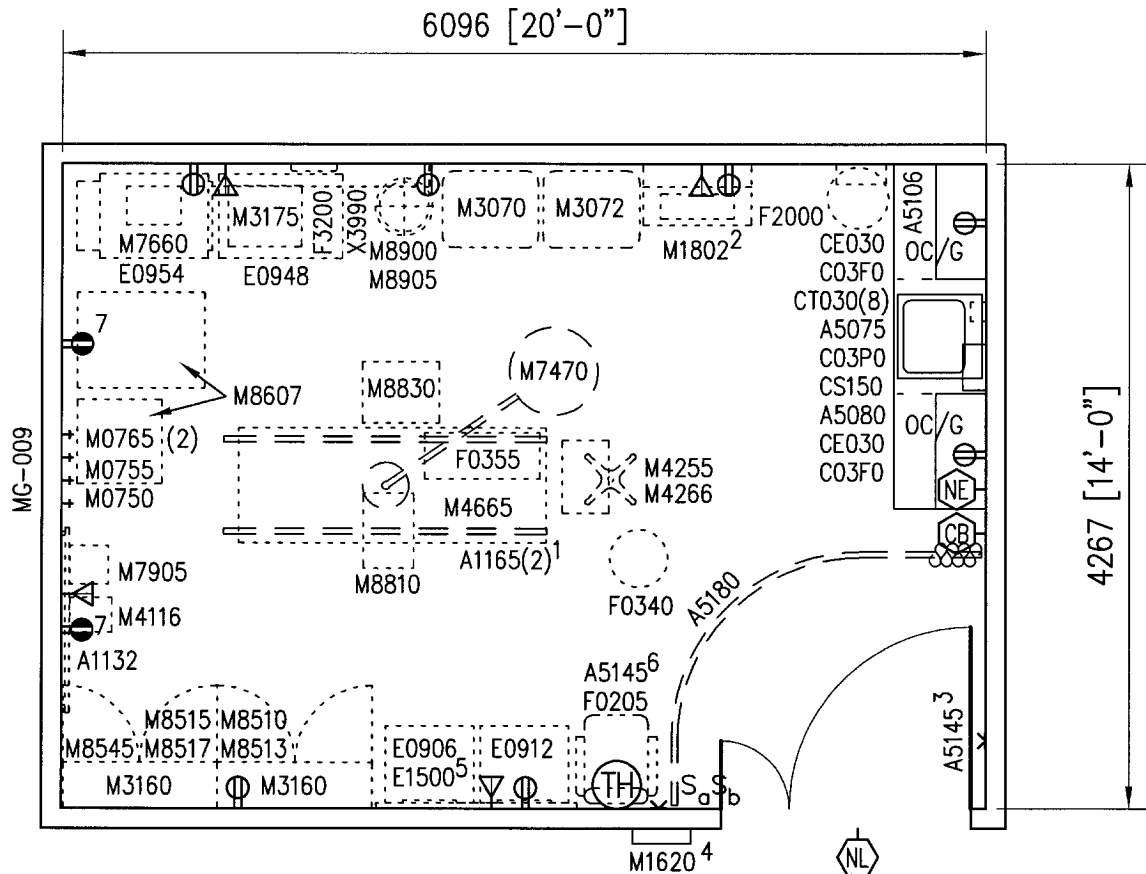

#### NOTES

1. PROVIDE X4880, MOBILE FLUOROSCOPIC C-ARM, IN VICINITY OF ROOM FOR NAVY PROJECTS ONLY.

#### KEY NOTES

1. PROVIDE A1165, IV TRACKS, FOR NAVY PROJECTS ONLY.
2. COORDINATE LOCATIONS OF ELECTRICAL AND DATA OUTLETS WITH M1802, WALL MOUNTED, FOLD DOWN COMPUTER.
3. MOUNT A5145, HOOK, AT 1372 [4'-6"] AFF. PROVIDE 2-A5145'S FOR NAVY PROJECTS ONLY. MOUNT 1-A5145 AT 1372 [4'-6"] AFF AND 1-A5145 AT 1524 [5'-0"] AFF.
4. COORDINATE LOCATION OF M1620, PATIENT CHART HOLDER, AND ROOM SIGNAGE.
5. LENGTH OF E1500, WALL MOUNTED RAIL, TO BE DETERMINED BY MANUFACTURER/INSTALLER.
6. MOUNT A5145, HOOK, AT 1219 [4'-0"] AFF.
7. UTILITIES AND DEVICES INDICATED ARE FOR A CLINIC OR HOSPITAL WITH AN ALTERNATE POWER SOURCE.

## TREATMENT ROOM, ENDOSCOPY EXAM (TREE1)

**3.15.04**

SPECIALTY MEDICAL CLINICS

FLOOR/EQUIPMENT PLAN


26.01 NSM (280 NSF)

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

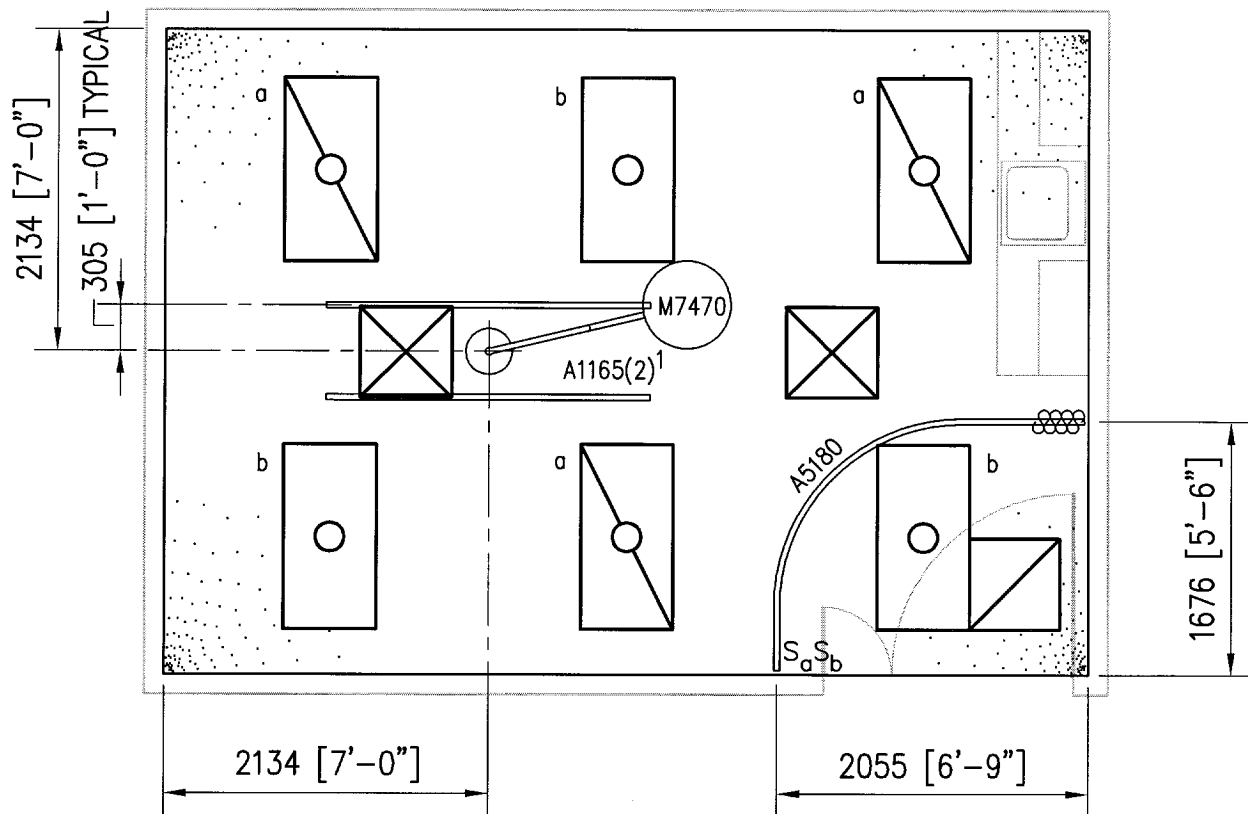

#### KEY NOTES

1. PROVIDE A1165, IV TRACKS, FOR NAVY PROJECTS ONLY.

## TREATMENT ROOM, ENDOSCOPY EXAM (TREE1)

**3.15.04**

SPECIALTY MEDICAL CLINICS

CEILING PLAN


DATE: 01/28/2002

0 1000 2000




1:50 SCALE (MILLIMETERS)



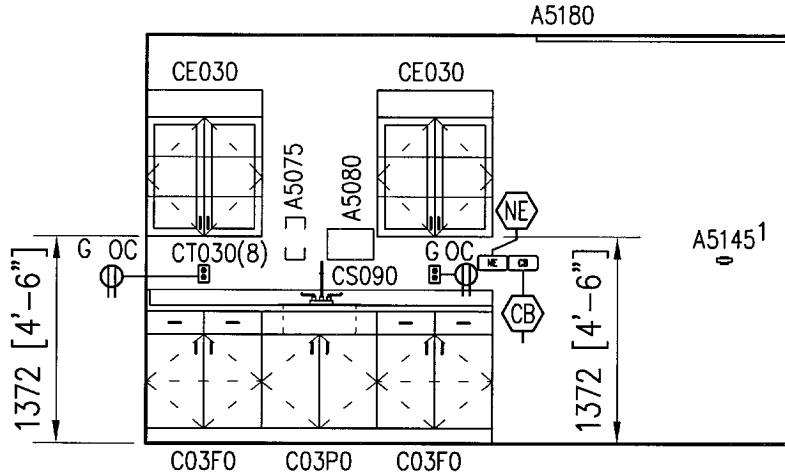

**EAST**

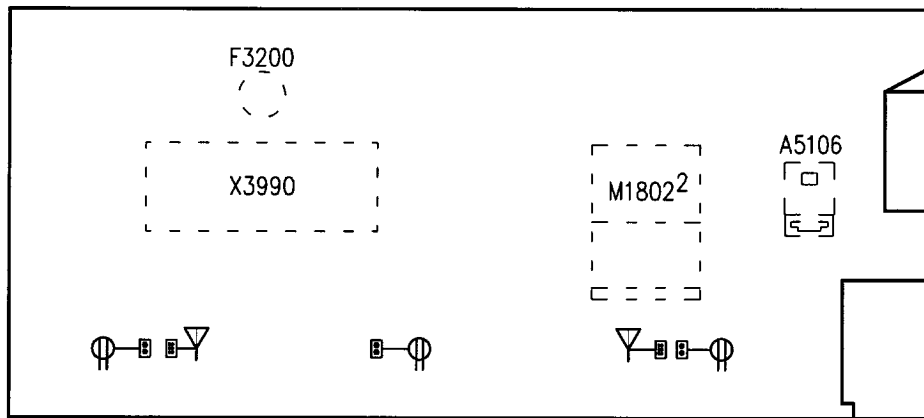

**NORTH**

**NOTES**

1. SEE TYPICAL EQUIPMENT MOUNTING HEIGHTS BEGINNING ON SHEET 7.01.01.

**KEY NOTES**

1. MOUNT A5145, HOOK, AT 1372 [4'-6"] AFF. PROVIDE 2-A5145'S FOR NAVY PROJECTS ONLY. MOUNT 1-A5145 AT 1372 [4'-6"] AFF AND 1-A5145 AT 1524 [5'-0"] AFF.
2. COORDINATE LOCATIONS OF ELECTRICAL AND DATA OUTLETS WITH M1802, WALL MOUNTED, FOLD DOWN COMPUTER.

**TREATMENT ROOM,  
ENDOSCOPY EXAM (TREE1)**

**3.15.04**

SPECIALTY MEDICAL CLINICS


ELEVATIONS

DATE: 01/28/2002


0 1000 2000



1:50 SCALE (MILLIMETERS)

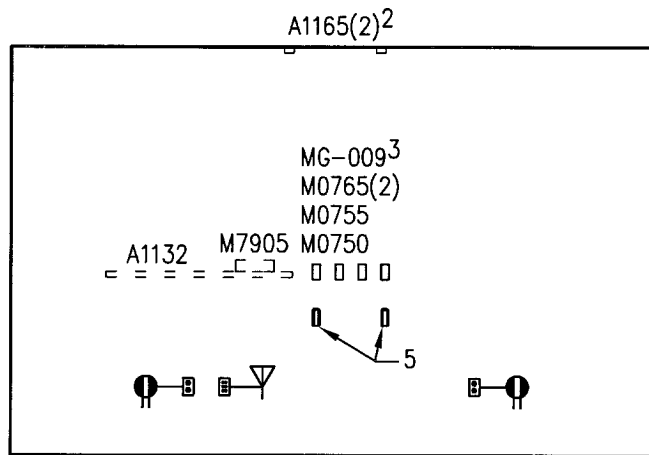

**WEST**

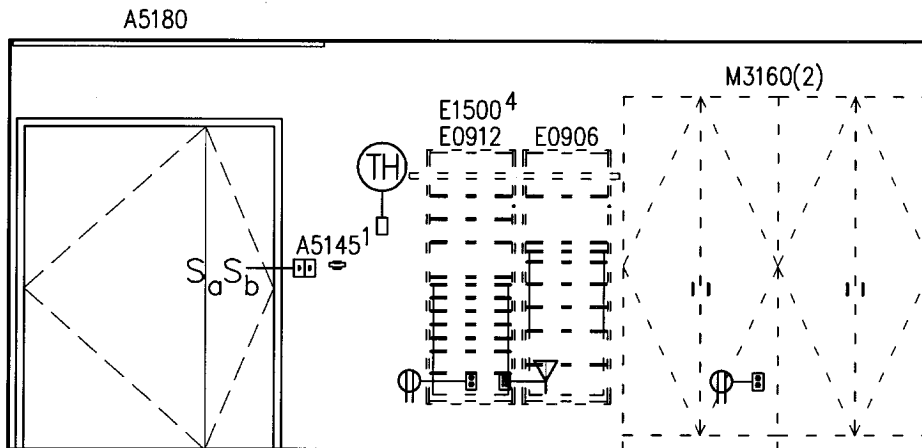

**SOUTH**

**NOTES**

1. SEE TYPICAL EQUIPMENT MOUNTING HEIGHTS BEGINNING ON SHEET 7.01.01.

**KEY NOTES**

1. MOUNT A5145, HOOK, AT 1219 [4'-0"] AFF.
2. PROVIDE A1165, IV TRACKS, FOR NAVY PROJECTS ONLY.
3. PROVIDE BOTTLED GAS IF BUILDING PIPED MCA AND OX ARE NOT AVAILABLE.
4. LENGTH OF E1500, WALL MOUNTED RAIL, TO BE DETERMINED BY MANUFACTURER/INSTALLER.
5. VACUUM SLIDE.

**TREATMENT ROOM,  
ENDOSCOPY EXAM (TREE1)**

**3.15.04**

SPECIALTY MEDICAL CLINICS

ELEVATIONS

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

# Medical Equipment Room Guide Plates

## 3.15 Specialty Medical Clinics

**Plate Number** 3.15.04


**Room Name** TREATMENT ROOM, ENDOSCOPY EXAM

**SEPS #** TREE1

JSN	Nomenclature	Qty	Army	Navy	AF	Utilities					
			Cat	Cat	Cat	1	2	3	4	5	6
A1132	Rail, Accessory Mounting, Length As Required	1	C	C	C	.	.	.	.	.	.
A1165	Track, IV, Ceiling Mounted, 7 Foot	2	A	A	A	.	.	.	.	.	J
A5075	Dispenser, Soap, Disposable	1	C	C	C	.	.	.	.	.	.
A5080	Dispenser, Paper Towel, SS, Surface Mounted	1	A	A	A	.	.	.	.	.	J
A5106	Waste Disposal Unit, Sharps w/Glove Dispenser	1	C	A	C	.	.	.	.	.	.
A5145	Hook, Garment, Double, SS, Surface Mounted	2	A	A	A	.	.	.	.	.	J
A5180	Track, Cubicle, Surface Mounted, With Curtain	1	A	A	A	.	.	.	.	.	J
C03F0	Cabinet, U/C/B, 1 Shelf, 2 Half DR, 2 DO, 36x30x22	2	A	A	A	.	.	.	.	.	.
C03P0	Cabinet, Sink, U/C/B, 2 Door, 36x30x22	1	A	A	A	.	.	.	.	.	.
CE030	Cabinet, W/H, 2 SH, 2 GDO, Sloping Top, 38x30x13	2	A	A	A	.	.	.	.	.	J
CS150	Sink, SS, Single Compartment, 10x19x16 ID	1	A	A	A	D	.	.	.	.	.
CT030	Countertop, High Pressure Laminate	8	A	A	A	.	.	.	.	.	.
E0906	Locker, Supply, General, Wall Mtd, 23"W x 20"D	1	C	C	C	.	.	.	.	.	.
E0912	Locker, Supply, Med Surg, Wall Mtd, 23"W x 20"D	1	C	C	C	.	.	.	.	.	.
E0948	Cart, General Storage, Mobile, 42"H x 32"W x 22"D	1	C	C	C	.	.	.	.	.	.
E0954	Cart, Emergency, Mobile, 66"H x 52"W x 22"D	1	C	C	C	.	.	.	.	.	.
E1500	Rail, MOD, W/MNTD, HX144XD	1	C	C	C	.	.	.	.	.	J
F0205	Chair, Side With Arms	1	C	C	C	.	.	.	.	.	.
F0340	Stool, Self Adjusting	1	C	C	C	.	.	.	.	.	.
F0355	Footstool, Straight	1	C	C	C	.	.	.	.	.	.
F2000	Basket, Wastepaper, Round, Metal, 18 H x 16 Dia.	1	C	C	C	.	.	.	.	.	.
F3200	Clock, Battery, 12" Diameter	1	C	C	C	.	.	.	.	.	.
M0750	Flowmeter, Air, Connect w/50 PSI Supply	1	C	C	C	.	.	K	.	.	.
M0755	Flowmeter, Oxygen, Low Flow	1	C	C	C	.	.	A	.	.	.
M0765	Regulator, Vacuum	2	C	C	C	.	.	B	.	.	.
M1620	Holder, Chart, Patient, Wall or Door Mounted	1	A	A	A	.	.	.	.	.	.
M1802	Wall Mounted, Fold Down Computer	1	C	C	C	.	A	.	.	.	J
M3070	Hamper, Linen, Mobile, w/Lid	1	C	C	C	.	.	.	.	.	.
M3072	Frame, Infectious Waste Bag w/Lid	1	C	C	C	.	.	.	.	.	.
M3160	Cabinet, Storage, Hanging Scopes	2	C	C	C	.	.	.	.	.	.
M3175	Electrosurgical Unit, Dual Output	1	C	C	C	.	A	.	.	.	.
M4116	Monitor, Vital Signs	1	C	C	C	.	A	.	.	.	.
M4255	Stand, IV, Adjustable	1	C	C	C	.	.	.	.	.	.
M4266	Pump, Volumetric, Infusion, Multiple Line	1	C	C	C	.	A	.	.	.	.
M4665	Stretcher, Recovery, Surgical	1	C	C	C	.	.	.	.	.	.
M7470	Light, Surgical, Ceiling, Single, Small	1	A	A	A	.	F	.	.	.	J
M7660	Defibrillator/Monitor/Recorder, Portable	1	C	C	C	.	A	.	.	.	.
M7905	Oximeter, Pulse	1	C	C	C	.	A	.	.	.	.
M8510	Colonoscope, Fiberoptic, Therapeutic	1	C	C	C	.	.	.	.	.	.

## Medical Equipment Room Guide Plates

M8513	Duodenoscope, Fiberoptic, Therapeutic	1	C	C	C	.	.	.	.	.	.
M8515	Gastroscope, Fiberoptic, Therapeutic	1	C	C	C	.	.	.	.	.	.
M8517	Enteroscope, Fiberoptic, Therapeutic	1	C	C	C	.	.	.	.	.	.
M8545	Sigmoidoscope, Fiberoptic, Therapeutic	1	C	C	C	.	.	.	.	.	.
M8607	Endoscopy Cart, Fiberoptic, w/Video, Info. Mgmt.	1	C	C	C	.	A	.	E	.	.
M8810	Stand, Mayo	1	C	C	C	.	.	.	.	.	.
M8830	Table, Instrument/Dressing, Mobile, 34x20x16	1	C	C	C	.	.	.	.	.	.
M8900	Carriage, Pail, CRS, Without Pail	1	C	C	C	.	.	.	.	.	.
M8905	Pail, Utility, CRS, Without Carriage	1	C	C	C	.	.	.	.	.	.
X3990	Illuminator, Film, 4 Panels, Wall Mounted	1	C	C	C	.	A	.	.	.	.

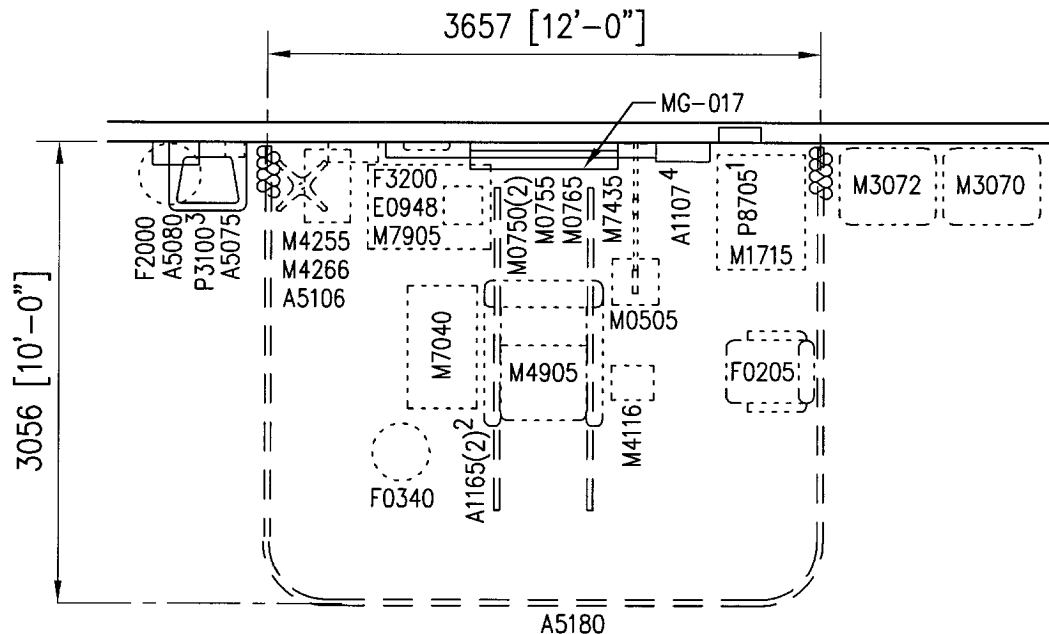

#### NOTES

1. AN E0954, EMERGENCY CART, AND A M7660, DEFIBRILLATOR, SHALL BE IN THE VICINITY OF THE DIALYSIS STATIONS.
2. M7850, CENTRAL PHYSIOLOGICAL MONITOR, AND M7845, BEDSIDE PHYSIOLOGICAL MONITOR, MAY BE ADDED. VERIFY WITH USER.
3. A M4030, SIT/STAND SCALE, AND A M3150, AUTOMATIC MEDICATION DISTRIBUTION SYSTEM, SHALL BE IN THE VICINITY OF THE DIALYSIS STATIONS.
4. TREATED WATER REQUIREMENTS MAY BE SATISFIED WITH EITHER A CENTRAL SYSTEM OR POINT-OF-USE SYSTEM. SYSTEM SELECTION SHALL BE COORDINATED WITH THE USER AND THE DESIGN AGENT.

#### KEY NOTES

1. P8705, DIALYSIS BOX, WITH COLD WATER AND DRAIN.
2. PROVIDE A1165, IV TRACKS, FOR NAVY PROJECTS ONLY.
3. P3200, ELECTRONIC LAVATORY, MAY BE USED INSTEAD OF P3100, LAVATORY. FOR THE P3200 PROVIDE EITHER BATTERY POWER OR HARD WIRED CONNECTION. PROVIDE ONE P3100 AND ACCESSORIES FOR EVERY 2 DIALYSIS STATIONS.
4. SEE ELEVATION OF A1107, UTILITY RAIL SYSTEM, FOR UTILITIES.

## RENAL DIALYSIS STATION, CHAIR (RDC01) 3.15.05

SPECIALTY MEDICAL CLINICS


FLOOR/EQUIPMENT PLAN 11.15 NSM (120 NSF)

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

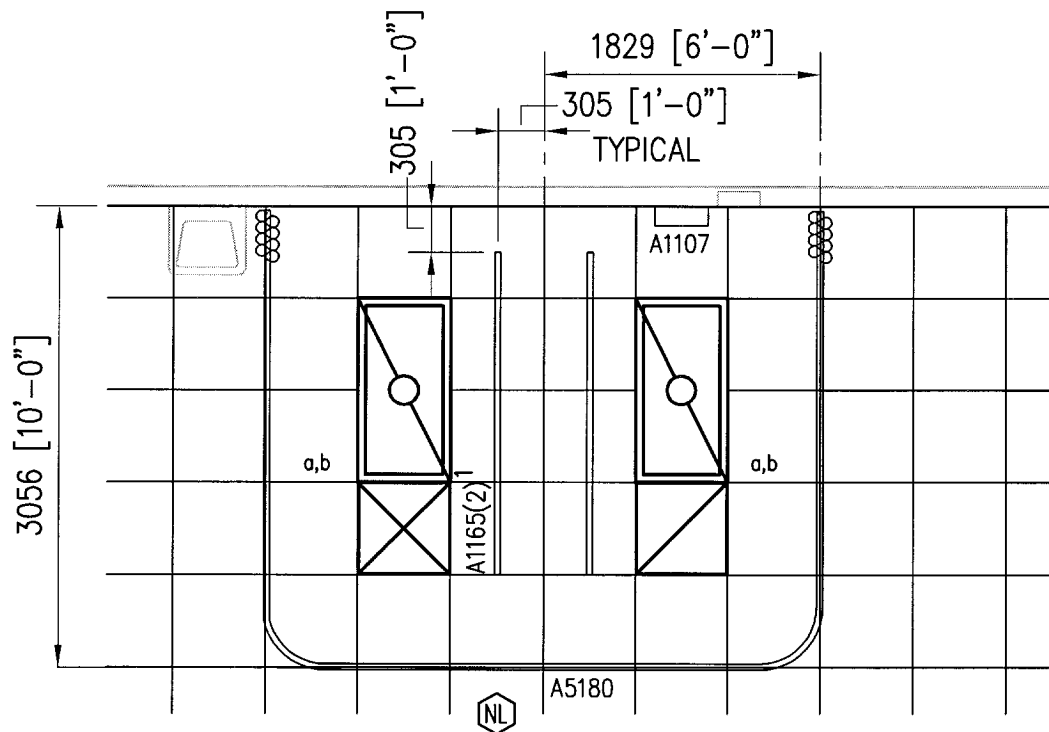

#### KEY NOTES

1. PROVIDE A1165, IV TRACKS, FOR NAVY PROJECTS ONLY.

## RENAL DIALYSIS STATION, CHAIR (RDC01) 3.15.05

SPECIALTY MEDICAL CLINICS


CEILING PLAN

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

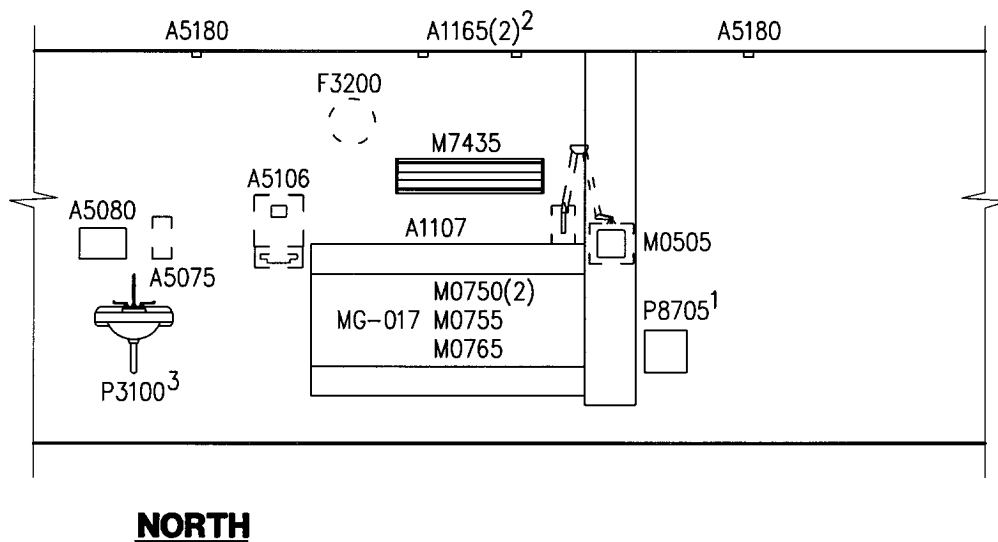

**NOTES**

1. SEE TYPICAL EQUIPMENT MOUNTING HEIGHTS BEGINNING ON SHEET 7.01.01.

**KEY NOTES**

1. P8705, DIALYSIS BOX, WITH COLD WATER AND DRAIN.
2. PROVIDE A1165, IV TRACKS, FOR NAVY PROJECTS ONLY.
3. P3200, ELECTRONIC LAVATORY, MAY BE USED INSTEAD OF P3100, LAVATORY. FOR THE P3200 PROVIDE EITHER BATTERY POWER OR HARD WIRED CONNECTION. PROVIDE ONE P3100 AND ACCESSORIES FOR EVERY 2 DIALYSIS STATIONS.

**RENAL DIALYSIS STATION, CHAIR (RDC01) 3.15.05**

SPECIALTY MEDICAL CLINICS

ELEVATIONS

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)






## Medical Equipment Room Guide Plates

### 3.15 Specialty Medical Clinics

**Plate Number** 3.15.05  
**Room Name** RENAL DIALYSIS STATION, CHAIR  
**SEPS #** RDC01

JSN	Nomenclature	Qty	Army	Navy	AF	Utilities					
			Cat	Cat	Cat	1	2	3	4	5	6
A1107	Rail System , Utility, Gas and Electric	1	A	A	A	.	F	G	.	.	A
A1165	Track, IV, Ceiling Mounted, 7 Foot	2	A	A	A	.	.	.	.	.	J
A5075	Dispenser, Soap, Disposable	1	C	C	C	.	.	.	.	.	.
A5080	Dispenser, Paper Towel, SS, Surface Mounted	1	A	A	A	.	.	.	.	.	J
A5106	Waste Disposal Unit, Sharps w/Glove Dispenser	1	C	A	C	.	.	.	.	.	.
A5180	Track, Cubicle, Surface Mounted, With Curtain	1	A	A	A	.	.	.	.	.	J
E0948	Cart, General Storage, Mobile, 42"H x 32"W x 22"D	1	C	C	C	.	.	.	.	.	.
F0205	Chair, Side With Arms	1	C	C	C	.	.	.	.	.	.
F0340	Stool, Self Adjusting	1	C	C	C	.	.	.	.	.	.
F2000	Basket, Wastepaper, Round, Metal, 18 H x 16 Dia.	1	C	C	C	.	.	.	.	.	.
F3200	Clock, Battery, 12" Diameter	1	C	C	C	.	.	.	.	.	.
M0505	Television, Color, Bedside	1	C	C	C	.	F	.	.	.	.
M0750	Flowmeter, Air, Connect w/50 PSI Supply	2	C	C	C	.	.	K	.	.	.
M0755	Flowmeter, Oxygen, Low Flow	1	C	C	C	.	.	A	.	.	.
M0765	Regulator, Vacuum	1	C	C	C	.	.	B	.	.	.
M1715	Hemodialysis Unit, Hollow Fiber or Plate	1	C	C	C	E	A	.	.	.	.
M3070	Hamper, Linen, Mobile, w/Lid	1	C	C	C	.	.	.	.	.	.
M3072	Frame, Infectious Waste Bag w/Lid	1	C	C	C	.	.	.	.	.	.
M4116	Monitor, Vital Signs	1	C	C	C	.	A	.	.	.	.
M4255	Stand, IV, Adjustable	1	C	C	C	.	.	.	.	.	.
M4266	Pump, Volumetric, Infusion, Multiple Line	1	C	C	C	.	A	.	.	.	.
M4905	Chair, Dialysis	1	C	C	C	.	.	.	.	.	.
M7040	Table, Overbed	1	C	C	C	.	.	.	.	.	.
M7435	Light, Overbed, Direct And Indirect	1	A	A	A	.	F	.	.	.	J
M7905	Oximeter, Pulse	1	C	C	C	.	A	.	.	.	.
P3100	Lavatory, Vitreous China, Slab Type	1	A	A	A	D	.	.	.	.	.
P8705	Dialysis Box, Recessed, Single Valve	1	A	A	A	E	.	.	.	.	.

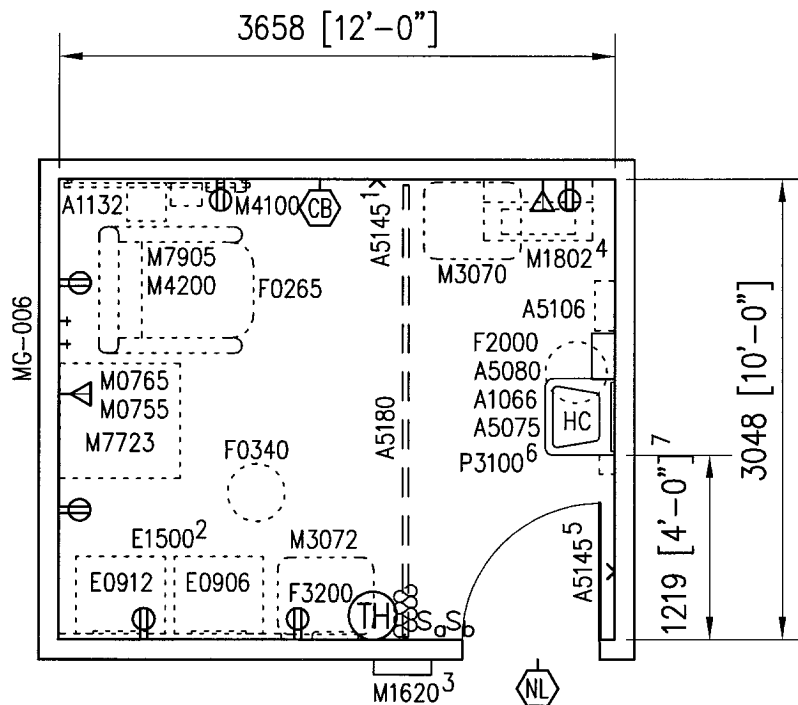

#### KEY NOTES

1. MOUNT A5145, HOOK, AT 1372 [4'-6"] AFF.
2. LENGTH OF E1500, WALL MOUNTED RAIL, TO BE DETERMINED BY MANUFACTURER/INSTALLER.
3. COORDINATE LOCATION OF M1620, PATIENT CHART HOLDER, AND ROOM SIGNAGE.
4. COORDINATE LOCATION OF ELECTRICAL AND DATA OUTLETS WITH M1802, WALL MOUNTED, FOLD DOWN COMPUTER.
5. MOUNT A5145, HOOK, AT 1219 [4'-0"] AFF. PROVIDE 2-A5145'S FOR NAVY PROJECTS ONLY. MOUNT 1-A5145 AT 1219 [4'-0"] AFF AND 1-A5145 AT 1524 [5'-0"] AFF.
6. P3200, ELECTRONIC LAVATORY, MAY BE USED INSTEAD OF P3100, LAVATORY. FOR THE P3200 PROVIDE EITHER BATTERY POWER OR HARD WIRED CONNECTION.
7. INCREASE THIS DIMENSION AS NEEDED TO COMPLY WITH ADA IF THE DOOR HAS A CLOSER OR IF 610 [2'-0"] CLEARANCE IS NOT AVAILABLE AT THE STRIKE SIDE OF THE DOOR.

## OUTPATIENT EEG TESTING AREA, 1 STATION (OPEE1) 3.15.06

SPECIALTY MEDICAL CLINICS


FLOOR/EQUIPMENT PLAN 11.15 NSM (120 NSF)

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

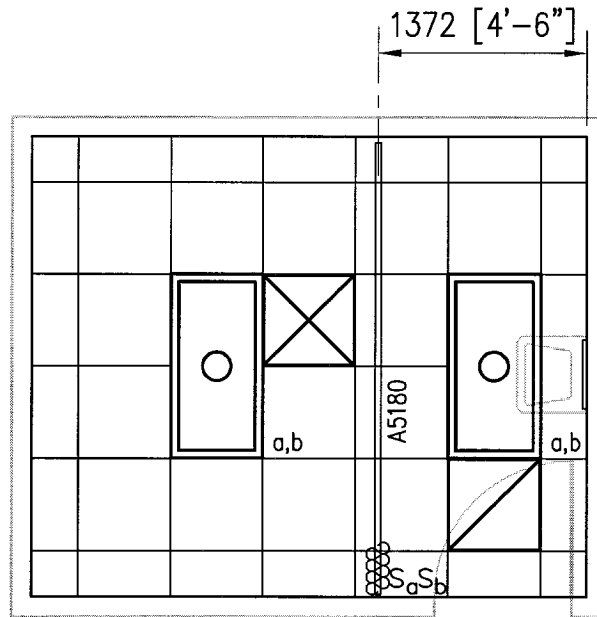

# **OUTPATIENT EEG TESTING AREA, 1 STATION (OPEE1) 3.15.06**

SPECIALTY MEDICAL CLINICS


CEILING PLAN

DATE: 01/28/2002


0 1000 2000




1:50 SCALE (MILLIMETERS)

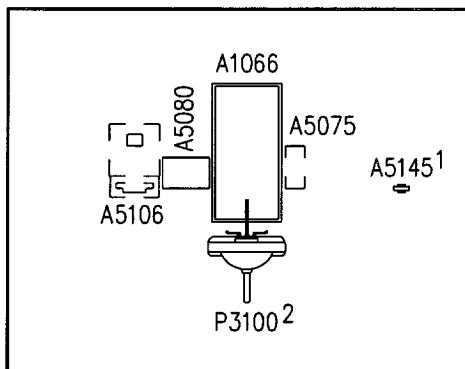


**EAST**

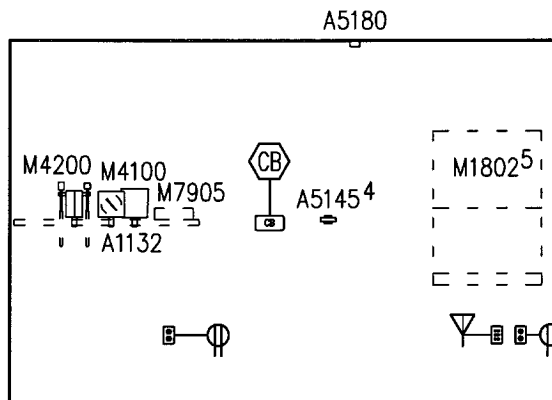

**NORTH**

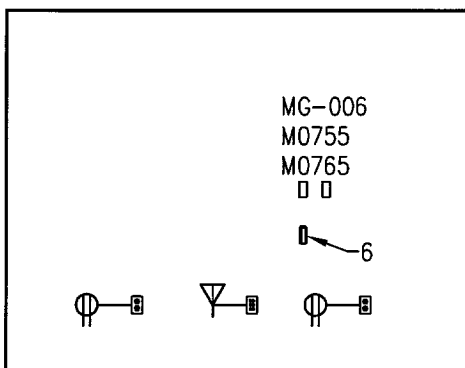

**WEST**

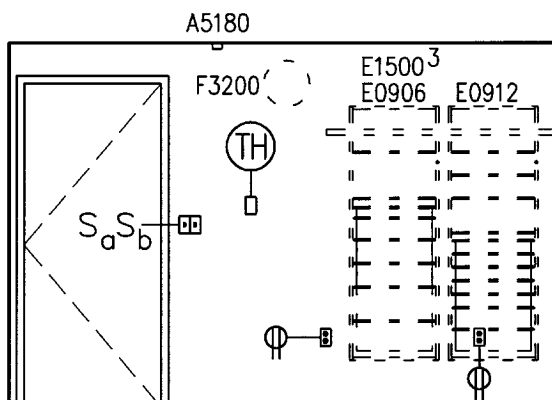

**SOUTH**

**NOTES**

1. SEE TYPICAL EQUIPMENT MOUNTING HEIGHTS BEGINNING ON SHEET 7.01.01.

**KEY NOTES**

1. MOUNT A5145, HOOK, AT 1219 [4'-0"] AFF. PROVIDE 2-A5145'S FOR NAVY PROJECTS ONLY. MOUNT 1-A5145 AT 1219 [4'-0"] AFF AND 1-A5145 AT 1524 [5'-0"] AFF.
2. P3200, ELECTRONIC LAVATORY, MAY BE USED INSTEAD OF P3100, LAVATORY. FOR THE P3200 PROVIDE EITHER BATTERY POWER OR HARD WIRED CONNECTION.
3. LENGTH OF E1500, WALL MOUNTED RAIL, TO BE DETERMINED BY MANUFACTURER/INSTALLER.
4. MOUNT A5145, HOOK, AT 1372 [4'-6"] AFF.
5. COORDINATE LOCATIONS OF ELECTRICAL AND DATA OUTLETS WITH M1802, WALL MOUNTED, FOLD DOWN COMPUTER.
6. VACUUM SLIDE.

# **OUTPATIENT EEG TESTING AREA, 1 STATION (OPEE1)**

**3.15.06**

SPECIALTY MEDICAL CLINICS

ELEVATIONS

DATE: 01/28/2002

0 1000 2000



1:50 SCALE (MILLIMETERS)

## Medical Equipment Room Guide Plates

### 3.15 Specialty Medical Clinics

**Plate Number** 3.15.06

**Room Name** OUTPATIENT EEG TESTING AREA, 1 STATION

**SEPS #** OPEE1

JSN	Nomenclature	Qty	Army	Navy	AF	Utilities					
			Cat	Cat	Cat	1	2	3	4	5	6
A1066	Mirror, Float Glass, With SS Frame, 36X18	1	A	A	A	.	.	.	.	.	J
A1132	Rail, Accessory Mounting, Length As Required	1	C	C	C	.	.	.	.	.	.
A5075	Dispenser, Soap, Disposable	1	C	C	C	.	.	.	.	.	.
A5080	Dispenser, Paper Towel, SS, Surface Mounted	1	A	A	A	.	.	.	.	.	J
A5106	Waste Disposal Unit, Sharps w/Glove Dispenser	1	C	A	C	.	.	.	.	.	.
A5145	Hook, Garment, Double, SS, Surface Mounted	2	A	A	A	.	.	.	.	.	J
A5180	Track, Cubicle, Surface Mounted, With Curtain	1	A	A	A	.	.	.	.	.	J
E0906	Locker, Supply, General, Wall Mtd, 23"W x 20"D	1	C	C	C	.	.	.	.	.	.
E0912	Locker, Supply, Med Surg, Wall Mtd, 23"W x 20"D	1	C	C	C	.	.	.	.	.	.
E1500	Rail, MOD, W/MNTD, HX144XD	1	C	C	C	.	.	.	.	.	J
F0265	Chair, Recliner	1	C	C	C	.	.	.	.	.	.
F0340	Stool, Self Adjusting	1	C	C	C	.	.	.	.	.	.
F2000	Basket, Wastepaper, Round, Metal, 18 H x 16 Dia.	1	C	C	C	.	.	.	.	.	.
F3200	Clock, Battery, 12" Diameter	1	C	C	C	.	.	.	.	.	.
M0755	Flowmeter, Oxygen, Low Flow	1	C	C	C	.	.	A	.	.	.
M0765	Regulator, Vacuum	1	C	C	C	.	.	B	.	.	.
M1620	Holder, Chart, Patient, Wall or Door Mounted	1	A	A	A	.	.	.	.	.	.
M1802	Wall Mounted, Fold Down Computer	1	C	C	C	.	A	.	.	.	J
M3070	Hamper, Linen, Mobile, w/Lid	1	C	C	C	.	.	.	.	.	.
M3072	Frame, Infectious Waste Bag w/Lid	1	C	C	C	.	.	.	.	.	.
M4100	Sphygmomanometer, Aneroid, Wall Mounted	1	C	C	C	.	.	.	.	.	.
M4200	Otoscope/Ophthalmoscope, Wall Mounted	1	C	C	C	.	A	.	.	.	.
M7723	Electroencephalograph, 32 Channel	1	C	C	C	.	A	.	.	.	.
M7905	Oximeter, Pulse	1	C	C	C	.	A	.	.	.	.
P3100	Lavatory, Vitreous China, Slab Type	1	A	A	A	D	.	.	.	.	.